

Rev 1

TABLE OF CONTENTS

TABLE OF CONTENTS.....	1
1- Overview	2
2- Procedure	2


Rev 1

1. Overview

This procedure describes how to see the SPT Result.

2. Procedure

1. Select [Utility] and Select [Report Manager]. Then, click [Start].

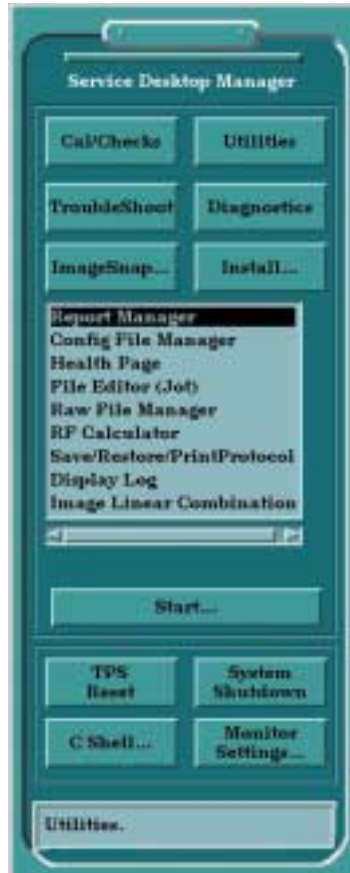

**SERVICE DESKTOP MANAGER
ILLUSTRATION 1**

Rev 1


2. Procedure(continued)

2. Enter “mrgo1” as password and click [Ok] button.



REPORT TOOL 1
ILLUSTRATION 2

3. Select “File / Open”. The following window comes up.
4. Select “SPT” from the File Type.
5. Select the file and click [OK] button.

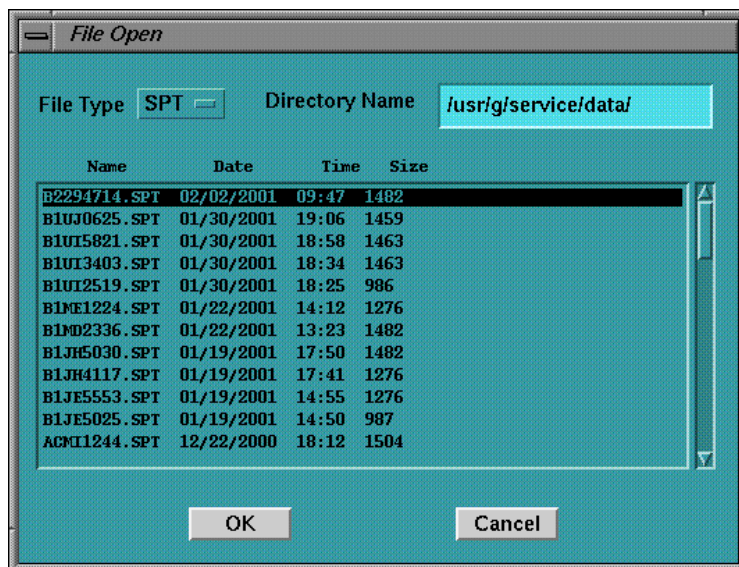


REPORT TOOL 2
ILLUSTRATION 3

Rev 1

2. Data Analysis (continued)


6. Header page is displayed. Click [Next] button.

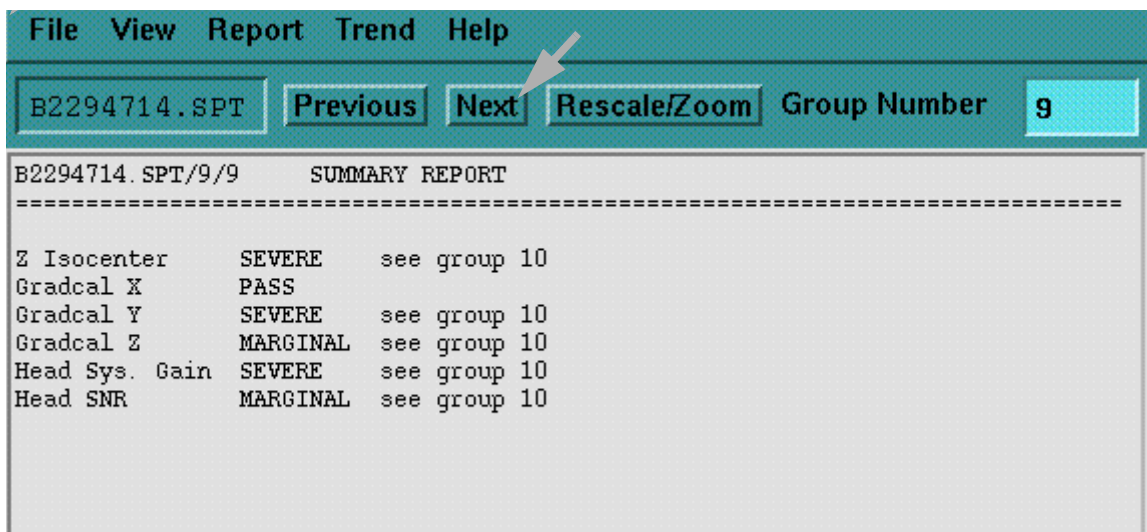

REPORT TOOL 3
ILLUSTRATION 4

7. The following screen comes up. Record "SNR". Click [Next] button.


Specification: SNR > 7.86

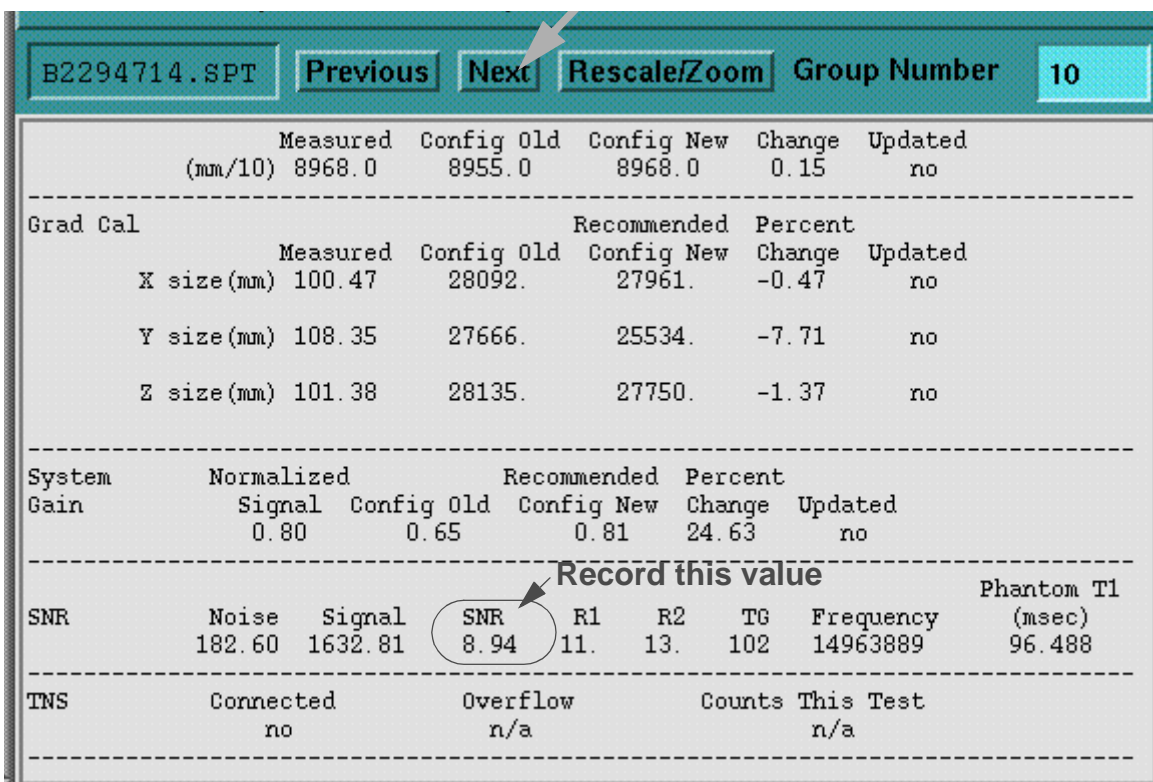

REPORT TOOL 4
ILLUSTRATION 5


Rev 1

8. The following screen comes up. Nothing to record. Click [Next] button.

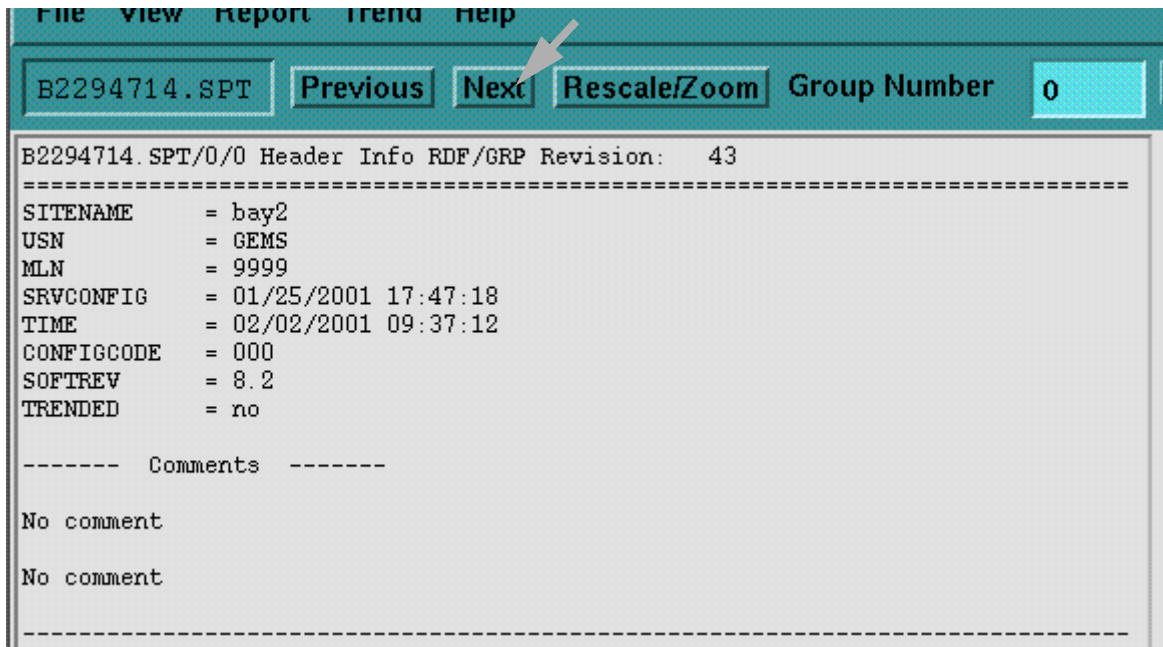

REPORT TOOL 4
ILLUSTRATION 6

9. The following screen comes up. Nothing to record.

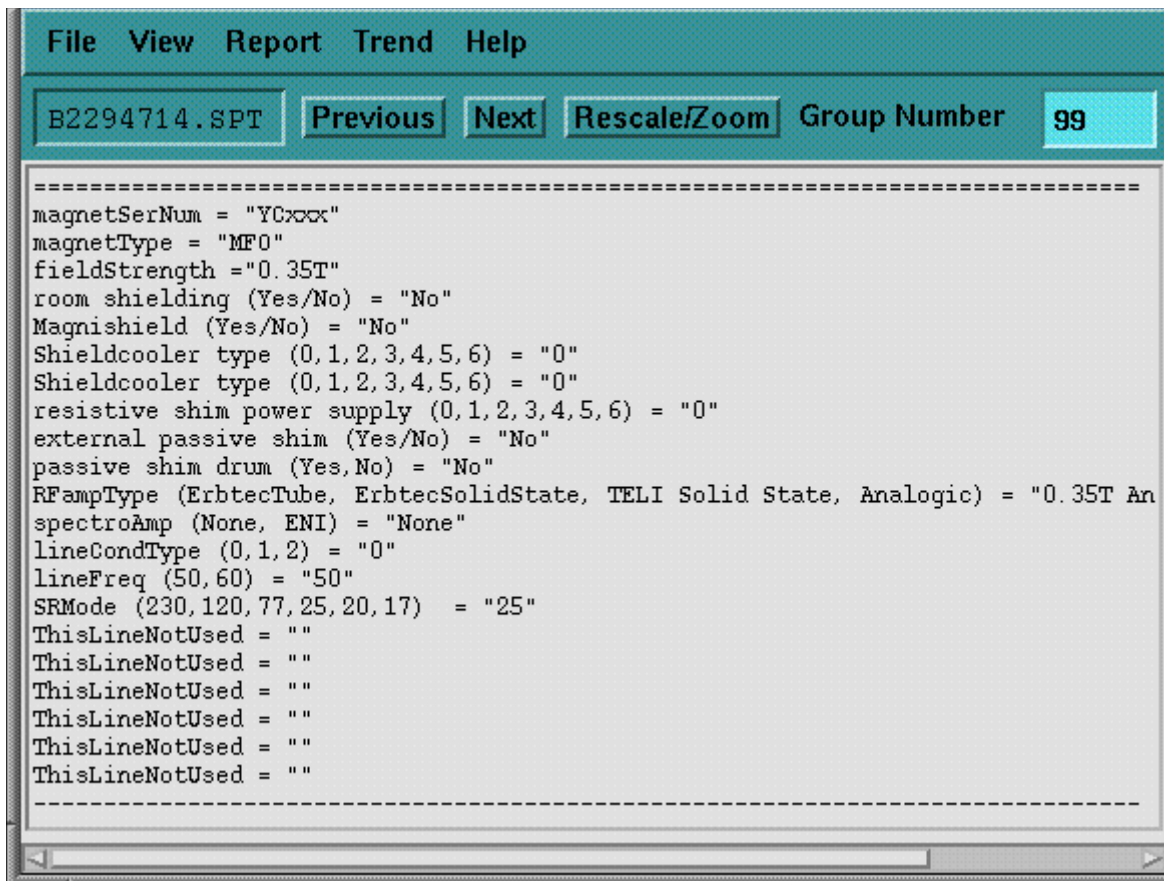

REPORT TOOL 5
ILLUSTRATION 7

Rev 1

Revision History

Rev	Date	Auther	Primary Reasons For Change
0	Feb 9, 2001	Y. Masumo	Initial Version
1	May 15, 2001	Y. Masumo	Added Specification